



BODYSTEEL & SILVER CO.,LTD
63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTABURI 11140 THAILAND

☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com

EO STERILIZATION Quality Control Report

Date : 21-Nov-2023

LOT : 0911

STERILIZATION DATE : 17-Nov-2023

EXPIRE DATE : 17-Nov-2028

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail




Yellow

Sample New Indicator



Pink

Quality Assurance By:

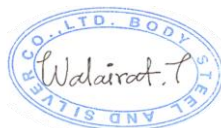

Date : **21-NOV-2023**

Welcome use sangiang medical equ Ethylene oxide sterilize		T °C		P kPa		H %		time	
Ster T:	50 °C	Start:	23.8	-0.3	89.7	16:50			
Ster H:	55 %	Water:	23.8	-0.2	74.6	16:51			
Ster time:	6h 0min	Vacuu:	23.8	-5.0	74.4	16:51			
Clean time:	5.0 min	Prhet:	35.0	-2.9	74.0	17:06			
Clean P:	-50 kPa	KHeat:	35.2	-2.0	74.4	17:16			
Clean time:	15 min	Vacuu:	34.0	-60.0	70.8	17:21			
		keepP:	40.7	-88.8	74.1	17:29			
		Inta:	39.7	-4.9	70.4	17:32			
		Vacuu:	39.3	-60.0	68.4	17:37			
		Humid:	39.3	-60.0	68.4	17:37			
		KHum:	38.8	-37.7	74.8	18:07			
		Dost:	38.1	-88.1	66.8	18:16			
		Rst:	38.2	-49.0	71.2	18:16			
		Diffu:	37.9	-24.6	70.8	18:26			
		Heat:	50.1	-14.3	74.0	19:02			
		Steri:	50.0	-14.3	74.1	19:02			
		Steri:	50.9	-11.5	74.1	19:12			
		Steri:	51.0	-10.5	74.9	19:22			
		Steri:	50.6	-10.1	74.0	19:32			
		Steri:	50.3	-9.6	74.3	19:42			
		Steri:	50.1	-9.1	74.7	19:52			
		Steri:	50.1	-8.8	74.9	20:02			
		Steri:	50.3	-8.4	74.7	20:12			
		Steri:	50.0	-8.0	74.2	20:22			
		Steri:	50.3	-7.7	74.5	20:32			
		Steri:	50.3	-7.4	74.5	20:42			
		Steri:	50.1	-7.1	74.9	20:52			
		Steri:	50.1	-6.8	74.5	21:02			
		Steri:	50.4	-6.5	74.2	21:12			
		Steri:	50.2	-6.2	74.3	21:22			
		Steri:	50.4	-5.9	74.0	21:32			
		Steri:	50.3	-5.6	74.9	21:42			
		Steri:	50.2	-5.4	74.3	21:52			
		Steri:	50.2	-5.1	74.9	22:02			
		Steri:	50.2	-4.8	74.5	22:12			
		Steri:	50.1	-4.4	74.4	22:22			
		Steri:	50.1	-4.1	74.1	22:32			
		Steri:	50.2	-3.7	74.9	22:42			
		Steri:	50.1	-3.4	74.7	22:52			
		Steri:	50.0	-3.0	74.3	23:02			
		Steri:	50.1	-2.7	74.9	23:12			
		Steri:	50.3	-2.4	74.6	23:22			
		Steri:	50.3	-2.0	74.3	23:32			
		Steri:	50.3	-1.8	74.3	23:42			
		Steri:	50.3	-1.6	74.4	23:52			
		Steri:	50.2	-1.5	74.5	00:02			
		Steri:	50.1	-1.3	74.3	00:12			
		Steri:	50.1	-1.1	74.3	00:22			
		Steri:	50.1	-1.0	74.7	00:32			
		Steri:	50.1	-1.0	74.3	00:42			
		Steri:	50.2	-0.9	74.1	00:52			
		Steri:	50.1	-0.8	74.8	01:02			
Unfinished Sterilization									
Steril Time: 8h21min									
Usecount: 00099									
Date: 2023-11-17									
Number: 50000450051									
Operator:									

Protocol: