



BODYSTEEL & SILVER CO.,LTD
63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTABURI 11140 THAILAND

☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com

EO STERILIZATION Quality Control Report

Date : 08-DEC-2023

LOT : 913

STERILIZATION DATE : 29-NOV-2023

EXPIRE DATE : 29-NOV-2028

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail



Yellow


Sample New Indicator



Pink

Quality Assurance By:

Date: 08-DEC-2023

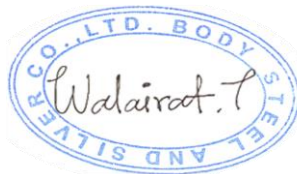

LOT NO. 913									
Welcome use sampling medical equ									
Ethylene oxide sterilize									
Ster T: 50 °C									
Ster H: 55 %									
SterTime: 6h 0min									
CleanTime: 5 min									
Clean P: -50 kPa									
CleanTime: 15 min									
T °C P kPa H % time									
Start:	22.7	-0.3	80.4	18:38					
Water:	22.9	-0.3	74.1	18:38					
Vacuum:	22.9	-5.1	74.9	18:38					
PHet:	35.0	-2.6	74.2	18:56					
KHeat:	35.2	-1.6	74.8	19:06					
Vacuum:	34.8	-60.1	62.9	19:11					
KeepP:	33.0	-69.0	74.0	19:19					
Inta:	38.5	-4.9	62.4	19:22					
Vacuum:	38.2	-60.1	80.1	19:27					
Humid:	38.2	-60.1	80.1	19:27					
KHum:	37.0	-57.8	74.5	19:57					
Dosi:	36.6	-58.1	61.9	20:06					
Bst:	36.8	-49.4	65.9	20:16					
Diffu:	37.7	-24.6	66.5	20:16					
Heat:	50.1	-18.0	74.8	20:42					
Steri:	50.0	-18.0	74.9	20:42					
Steri:	49.6	-15.7	74.1	20:52					
Steri:	51.5	-13.1	74.7	21:02					
Steri:	50.5	-10.8	74.7	21:12					
Steri:	50.7	-9.3	74.7	21:22					
Steri:	50.3	-8.6	74.3	21:32					
Steri:	50.2	-8.1	74.7	21:42					
Steri:	50.4	-7.5	74.9	21:52					
Steri:	50.1	-7.0	74.1	22:02					
Steri:	50.4	-6.4	74.1	22:12					
Steri:	50.2	-5.9	74.9	22:22					
Steri:	50.2	-5.4	74.7	22:32					
Steri:	50.2	-4.7	74.7	22:42					
Steri:	50.4	-4.8	74.1	22:52					
Steri:	50.3	-3.6	74.3	23:02					
Steri:	50.1	-2.9	74.9	23:12					
Steri:	50.0	-2.3	74.1	23:22					
Steri:	50.2	-1.8	74.9	23:32					
Steri:	50.2	-1.3	74.3	23:42					
Steri:	50.1	-1.0	74.9	23:52					
Steri:	50.4	-0.8	74.3	00:02					
Steri:	50.2	-0.6	74.9	00:12					
Steri:	50.2	-0.5	74.1	00:22					
Steri:	50.1	-0.5	74.5	00:32					
Steri:	50.2	-0.4	74.1	00:42					
Steri:	50.1	-0.3	74.1	00:52					
Steri:	50.2	-0.3	74.5	01:02					
Steri:	50.4	-0.3	74.5	01:12					
Steri:	50.2	-0.3	74.5	01:22					
Steri:	50.4	-0.3	74.1	01:32					
Steri:	50.2	-0.3	74.5	01:42					
Steri:	50.0	-0.4	74.7	01:52					
Steri:	50.1	-0.4	74.1	02:02					
Steri:	50.2	-0.4	74.7	02:12					
Steri:	50.0	-0.4	74.7	02:22					
Steri:	50.2	-0.4	74.7	02:32					
Steri:	50.3	-0.4	74.7	02:42					
Clean Process:	16.8	52.9	03:14						
PH:	50.6	-18.4	49.8	03:46					
2PH:	50.6	-18.4	49.8	03:46					
3PH:	50.5	-18.1	49.3	04:17					
4PH:	50.7	-18.2	48.3	04:49					
5PH:	50.7	-16.2	47.7	05:20					
Inta:	50.7	-9.5	45.6	06:21					
End:	50.7	-9.5	45.6	06:21					
Sterilization Completed									
Date: 2023-11-29 00:00									
Operator: 9000450051									

Protocol: