

**BODYSTEEL & SILVER CO.,LTD**63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTABURI 11140 THAILAND

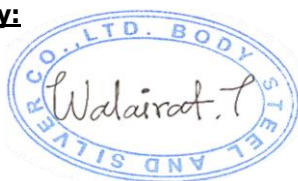
☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com**EO STERILIZATION
Quality Control Report****Date :15-DEC-2028****LOT : 915****STERILIZATION DATE : 9-DEC-2023****EXPIRE DATE : 9-DEC-2028****QC REPORT :****PASS as the test result has been shown as pictures below.****BIOLOGICAL INDICATORS :**

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

**Result
PASS****Sample Fail****Yellow****Sample New Indicator****Pink****Quality Assurance By:****Date: 15-DEC-2028**STERI 445 C DAT 0915
09/12/2023 16:52

STEP	TIME	TEMP(C)	PRES(KPa)	HUMID(CAR)
PH	16:52	28.8	00	84.5
PV	18:23	45.0	00	86.1
HU	18:36	44.8	-85	78.1
GS	18:46	44.9	-84	80.5
ST	18:48	45.6	-49	NONE
ST	19:18	45.2	-47	NONE
ST	19:48	44.6	-47	NONE
ST	20:18	45.1	-47	NONE
ST	20:48	45.1	-47	NONE
ST	21:18	44.7	-47	NONE
ST	21:48	45.5	-47	NONE
ST	22:18	44.8	-47	NONE
ST	22:48	44.8	-47	NONE
ST	23:18	45.4	-47	NONE
AE	23:48	44.7	-47	NONE
AE	00:47	45.3	-87	NONE
AE	01:47	45.1	-87	NONE
AE	02:48	45.3	-88	NONE
AE	03:48	44.9	-87	NONE
AE	04:48	45.4	-88	NONE
AE	05:47	44.9	-88	NONE
AE	06:47	45.3	-87	NONE
AE	07:46	44.9	-88	NONE
AE	08:44	45.4	-88	NONE
AE	09:43	44.9	-88	NONE
FI	09:48	45.5	-49	NONE

CHK BY :
APP BY :
RESULT :**Protocol:**