



BODYSTEEL & SILVER CO.,LTD
63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTABURI 11140 THAILAND

☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com

EO STERILIZATION Quality Control Report

Date : **12-JAN-2024**

LOT : 917

STERILIZATION DATE : 25-DEC-2023

EXPIRE DATE : 25-DEC-2028

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail



Yellow


Sample New Indicator



Pink

Quality Assurance By:

Date: 12-JAN-2024

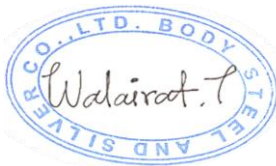

-STERI +45 C BAT :0917 25/12/2023 16:59		STEP	TIME	TEMP(C)	PRES(KPa)	HEAT(COOL)
PH	16:59	25.2	00	88.3		
PV	19:02	45.0	00	88.9		
HU	19:15	44.5	-85	89.7		
65	19:25	44.7	-84	83.9		
ST	19:27	45.6	-53	NONE		
ST	19:57	45.3	-52	NONE		
ST	20:27	44.7	-52	NONE		
ST	20:57	44.6	-52	NONE		
ST	21:27	45.3	-52	NONE		
ST	21:57	45.0	-51	NONE		
ST	22:27	44.4	-51	NONE		
ST	22:57	44.9	-51	NONE		
ST	23:27	45.4	-51	NONE		
ST	23:57	44.8	-51	NONE		
AE	00:27	44.6	-51	NONE		
AE	01:25	45.6	-88	NONE		
AE	02:23	44.8	-88	NONE		
AE	03:21	45.6	-88	NONE		
AE	04:19	44.7	-87	NONE		
AE	05:17	45.6	-88	NONE		
AE	06:15	44.8	-88	NONE		
AE	07:12	45.6	-88	NONE		
AE	08:09	44.8	-88	NONE		
AE	09:06	45.6	-87	NONE		
AE	10:02	44.8	-88	NONE		
FI	10:28	44.9	-88	NONE		
CHK BY :						
APP BY :						
RESULT :						

Protocol: