



BODYSTEEL & SILVER CO.,LTD
63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTABURI 11140 THAILAND

☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com

EO STERILIZATION Quality Control Report

Date : 17-May-2024

LOT : 0948

STERILIZATION DATE : 10-May-2024

EXPIRE DATE : 10-May-2029

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail




Yellow

Sample New Indicator



Pink

Quality Assurance By:

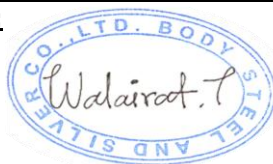

Date: 17-May-2024

lot 948		We know our long standing medical equip		
Phylone white autoclave				
Ster T:	50 °C			
Ster t:	55 min			
Ster line:	55 min			
Clean line:	5 min			
Clean p:	0.05			
Clean line:	15 min			
T	P	IF	time	
Ster1: 50.2	6.8	74.4	19:24	
Ster1: 50.3	6.2	74.4	19:34	
Ster1: 51.4	5.8	74.6	19:44	
Ster1: 50.6	5.3	74.2	19:54	
Ster1: 50.4	4.8	74.6	20:04	
Ster1: 50.1	4.3	74.4	20:14	
Ster1: 50.2	3.7	74.2	20:24	
Ster1: 50.0	3.1	74.4	20:34	
Ster1: 50.2	2.5	74.8	20:44	
Ster1: 50.3	2.0	74.4	20:54	
Ster1: 50.1	1.6	74.4	21:04	
Ster1: 50.1	1.2	74.4	21:14	
Ster1: 50.2	1.0	74.4	21:24	
Ster1: 50.0	0.8	74.2	21:34	
Ster1: 50.2	0.6	74.4	21:44	
Ster1: 50.2	0.5	74.2	21:54	
Ster1: 50.0	0.5	74.0	22:04	
Ster1: 50.3	0.4	74.8	22:14	
Ster1: 50.3	0.4	74.6	22:24	
Ster1: 50.2	0.1	74.8	22:34	
Ster1: 50.3	0.4	74.4	22:44	
Ster1: 50.2	0.4	74.2	22:54	
Ster1: 50.2	0.4	74.0	23:04	
Ster1: 50.1	0.4	74.6	23:14	
Ster1: 50.1	0.4	74.6	23:24	
Ster1: 50.1	0.4	74.4	23:34	
Ster1: 50.3	0.4	74.4	23:44	
Ster1: 50.3	0.4	74.4	23:54	
Ster1: 50.2	0.4	74.4	00:04	
Ster1: 50.2	0.4	74.0	00:17	
CLean Process				
1th :	51.0	17.5	74.3	00:30
2th :	50.6	17.4	74.0	01:05
3th :	50.5	17.3	74.4	01:32
4th :	50.4	17.4	74.1	01:59
5th :	50.4	17.4	74.0	02:26
Info :	50.4	10.4	74.3	02:27
End :	50.4	10.4	74.3	02:27
Sterilization Completed				
Ster1 time:	747min			
Use count:	0123			
Date:	2024-5-10			
Batch:	2024-5-10			
Operator:				

Protocol: