

EO STERILIZATION Quality Control Report

Date : 6-Sep-2024

LOT : 0984

STERILIZATION DATE : 3-Sep-2024

EXPIRE DATE : 3-Sep-2029

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail




Yellow

Sample New
Indicator

Pink

Quality Assurance By:

Date: 6-Sep-2024

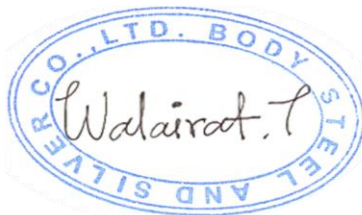

Lot 784

Welcome ase Sanitary medical ase
Ethylene oxide sterilize

Ster H: 55 °C
Ster L: 50 °C

Start time: 6h 00min
Clean time: 5 p/c
Clean P: 0 kPa
Clean time: 15 min

	T °C	P kPa	H %	time
Start:	46.2	-0.1	96.9	16:42
Water:	46.2	-0.0	74.7	16:43
Vacu:	46.0	-5.1	74.8	16:43
Priat:	46.0	-5.1	74.8	16:43
Kheat:	44.6	-2.2	74.2	16:53
Keap:	42.9	-68.8	74.6	16:58
Inta:	43.5	-4.9	74.9	17:10
Vacu:	42.0	-60.1	74.8	17:15
Humid:	42.0	-60.1	74.8	17:15
Humid:	40.6	-50.3	74.0	17:45
Dosi:	40.4	-57.5	74.0	17:53
Rst:	40.5	-48.8	74.8	17:54
Diffu:	40.1	-23.6	74.8	18:04
Heat:	50.1	-15.1	74.4	18:24
Ster:	50.1	-15.1	74.4	18:24
Ster:	49.7	-7.5	74.6	18:34
Ster:	50.0	-7.5	74.4	18:41
Ster:	50.6	-4.5	74.0	19:54
Ster:	50.6	-3.0	74.4	19:04
Ster:	50.4	-1.9	74.0	19:24
Ster:	50.5	-1.0	74.0	19:24
Ster:	50.2	-0.6	74.2	19:34
Ster:	50.2	-0.5	74.8	19:44
Ster:	50.2	-0.3	74.6	19:54
Ster:	50.1	-0.3	74.2	20:04
Ster:	50.1	-0.3	74.6	20:24
Ster:	50.2	-0.3	74.2	20:34
Ster:	50.0	-0.3	74.8	20:44
Ster:	50.0	-0.3	74.8	20:54
Ster:	50.5	-0.4	74.6	21:04
Ster:	50.2	-0.4	74.0	21:14
Ster:	50.2	-0.4	74.8	21:24
Ster:	50.2	-0.4	74.2	21:34
Ster:	50.3	-0.5	74.8	21:44
Ster:	50.0	-0.5	74.8	21:54
Ster:	50.2	-0.5	74.8	22:04
Ster:	50.2	-0.5	74.6	22:14
Ster:	50.4	-0.5	74.2	22:24
Ster:	50.2	-0.5	74.0	22:34
Ster:	50.2	-0.5	74.6	22:44
Ster:	50.2	-0.5	74.8	22:54
Ster:	50.2	-0.6	74.0	23:04
Ster:	50.3	-0.6	74.2	23:14
Ster:	50.1	-0.6	74.0	23:24
Ster:	50.1	-0.6	74.0	23:34
Ster:	50.4	-0.6	74.0	23:44
Ster:	50.3	-0.6	74.2	23:54
Ster:	50.3	-0.6	74.0	00:04
Ster:	50.2	-0.6	74.8	00:14
Ster:	50.1	-0.6	74.6	00:24
Clean	50.3	-11.3	74.7	00:50
1st:	50.2	-17.2	74.1	01:13
2nd:	50.3	-17.2	74.0	01:43
3rd:	50.2	-17.2	74.1	02:13
4th:	50.3	-17.1	74.3	02:37
5th:	50.3	-17.1	74.1	02:57
Inta:	50.4	-10.4	74.1	02:57
End:	50.4	-10.4	74.1	02:57

Sterilization Completed

Sterilization Monitor

Protocol: