



BODYSTEEL & SILVER CO.,LTD
63 Moo 5, Soi Wat Bang-Co Yaek Rong-Foam, Kanchanaphisek Road
Bang Ma-Nang, Bang Yai, NONTHABURI 11140 THAILAND

☎ 66 2-1505999 📠 66-2-1505990

sales@steelsilver.com

EO STERILIZATION Quality Control Report

Date : **25-APR-2025**

LOT : 1077

STERILIZATION DATE : 24-APR-2025

EXPIRE DATE : 24-APR-2029

QC REPORT :

PASS as the test result has been shown as pictures below.

BIOLOGICAL INDICATORS :

The EZTest indicator is placed in a suitable test pack which is representative of the load. The load is then processed as usual.

After sterilization the chemical indicator ring on the label changes from blue to brownish when exposed to Ethylene Oxide (EO gas). This distinguishes exposed from unexposed units (it does not indicate acceptable sterilization).

After the incubation period of 48 hours, the appearance of a yellow color indicates bacterial growth. No color change indicates adequate sterilization.

Result
PASS



Sample Fail



Yellow


Sample New Indicator



Pink

Quality Assurance By:

Date: 25-APR-2025

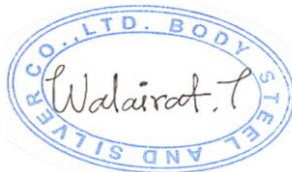

Welcome use sanjiang medical equipment ethylene oxide sterilize		T : P kPa H % time	
Lot 1099	Start: 24.1	-0.4	91.9 16:44
	Water: 24.2	-0.4	74.7 16:44
	Vacuum: 24.2	-5.2	74.8 16:45
	Pressure: 35.0	-2.7	74.7 17:00
	Heat: 35.0	-1.8	74.5 17:10
	Vacuum: 34.5	-60.0	74.6 17:15
	Temp: 38.9	-57.6	74.2 17:23
	Inta: 38.1	-4.9	74.8 17:26
	Vacuum: 38.4	-60.1	74.1 17:31
	Humid: 38.4	-60.1	74.3 17:31
	Kham: 37.6	-52.7	74.3 18:01
	Rest: 38.8	-58.2	74.3 18:08
	Diffus: 37.9	-74.0	74.8 18:19
	Heat: 50.1	-12.0	74.3 18:55
	Steri: 50.1	-8.9	74.9 19:05
	Steri: 50.1	-6.1	74.1 19:15
	Steri: 50.4	-3.9	74.3 19:25
	Steri: 50.2	-2.8	74.7 19:35
	Steri: 50.2	-1.9	74.5 19:45
	Steri: 50.2	-1.2	74.1 19:55
	Steri: 50.2	-0.8	74.3 20:05
	Steri: 50.0	-0.5	74.5 20:15
	Steri: 50.3	-0.3	74.3 20:25
	Steri: 50.2	-0.2	74.7 20:35
	Steri: 50.0	-0.2	74.7 20:45
	Steri: 50.5	-0.1	74.5 20:55
	Steri: 50.1	-0.1	74.9 21:05
	Steri: 50.2	-0.3	74.1 21:15
	Steri: 50.2	-0.2	74.1 21:25
	Steri: 50.2	-0.2	74.1 21:35
	Steri: 50.1	-0.2	74.7 21:45
	Steri: 50.1	-0.2	74.7 21:55
	Steri: 50.3	-0.2	74.5 22:05
	Steri: 50.2	-0.3	74.5 22:15
	Steri: 50.1	-0.3	74.3 22:25
	Steri: 50.2	-0.3	74.9 22:35
	Steri: 50.2	-0.3	74.9 22:45
	Steri: 50.1	-0.3	74.1 22:55
	Steri: 50.3	-0.3	74.5 23:05
	Steri: 50.3	-0.3	74.1 23:15
	Steri: 50.2	-0.3	74.3 23:25
	Steri: 50.0	-0.4	74.9 23:35
	Steri: 50.2	-0.4	74.7 23:45
	Steri: 50.1	-0.5	74.1 23:55
	Steri: 50.2	-0.4	74.1 00:05
	Steri: 50.0	-0.5	74.5 00:15
	Steri: 50.1	-0.5	74.1 00:25
	Steri: 50.2	-0.5	74.7 00:35
	Steri: 50.0	-0.5	74.1 00:45
	Steri: 50.0	-0.5	74.5 00:55
Clean Process			
	1th: 50.6	-15.1	74.8 01:21
	2th: 50.6	-15.3	74.2 01:47
	3th: 50.5	-15.2	74.2 02:14
	4th: 50.4	-15.1	74.2 02:44
	5th: 50.4	-15.2	74.4 03:07
	6th: 50.4	-15.2	74.1 03:34
	Inta: 50.4	-8.4	74.1 03:34
	End: 50.4	-8.4	74.1 03:34
Sterilization Completed			
	Steri Time: 10h50min		
	Usecount: 00183		
	Date: 2025-4-24		
	Number: 5000045051		
	Operator:		

Protocol: